

OLYMPUS

Your Vision, Our Future

OLYMPUS
EndoTherapy
BEYOND THE SCOPE

**Maximum force, minimum hassle— today's most reliable
and effective injection performance**

Single Use Injector

Injector
Force
Max™



Thanks to improved needle extrusion capability, enhanced click response, and excellent liquid feeding, the InjectorForce Max™ maximizes endoscopic injection performance throughout the gastrointestinal tract


Improved Extrusion Performance



Enanced piercing capability with improved extrusion performance delivers a smooth, reliable, effective injection every time

In addition to featuring an extraordinarily sharp needle tip, the InjectorForce Max™ offers much more effective extrusion performance over previous models. The combination of the sharp tip and improved extrusion performance provide just the right injection you need.

Improved Liquid Feeding

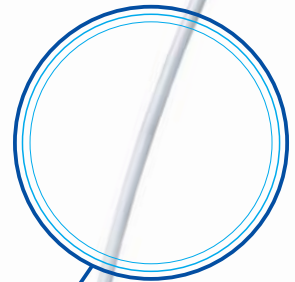



Improved Click Response


New ergonomic handle design and enhanced click response allow you to manipulate the needle tip more easily and accurately

Featuring an all-new ergonomic handle design and enhanced click response when you try to extrude the needle, the InjectorForce Max™ enables you to work with comfort and confidence.




A larger lumen ensures that the injection solution flows smoothly and quickly through the needle lumen

The InjectorForce Max™'s lumen is larger than the previous Injector Force's lumen to enable fluid to flow more smoothly and quickly, helping you to perform injections more efficiently and effectively than ever before.

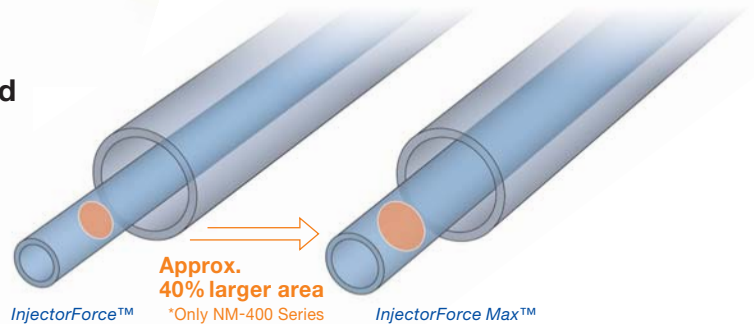

Wide lineup of versatile needles supports applications throughout the gastrointestinal tract

The InjectorForce Max™ wide-ranging lineup includes models with various needle diameters that can be passed through various scope channel diameters. The new design also includes an extensive selection of needle lengths and sheath lengths for use throughout the gastrointestinal tract from the esophagus and stomach to the small intestine and colon.

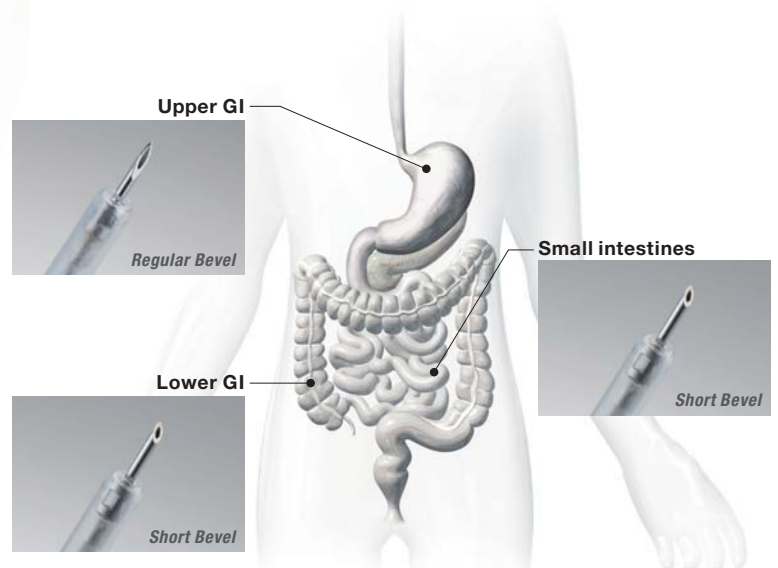

Specifications

	Model	Compatible Channel	Working Length	Bevel	Needle Length	Gauge
Upper GI	NM-400L-0421	2.8 mm	1650 mm	Regular	4 mm	21 G
	NM-400L-0423				4 mm	23 G
	NM-400L-0425				4 mm	25 G
	NM-400L-0523				5 mm	23 G
	NM-400L-0525				5 mm	25 G
	NM-400L-0621				6 mm	21 G
	NM-400L-0623				6 mm	23 G
	NM-400L-0625				6 mm	25 G
	NM-401L-0423	2.0 mm			4 mm	23 G
	NM-401L-0425				4 mm	25 G
	NM-401L-0523				5 mm	23 G
	NM-401L-0525				5 mm	25 G
	NM-401L-0623				6 mm	23 G
	NM-401L-0625				6 mm	25 G
Lower GI	NM-400U-0323	2.8 mm	2300 mm	Short	3 mm	23 G
	NM-400U-0423				4 mm	23 G
	NM-400U-0425				4 mm	25 G
	NM-400U-0523				5 mm	23 G
	NM-400U-0525				5 mm	25 G
	NM-400U-0623				6 mm	23 G
	NM-400U-0625				6 mm	25 G
For Small Intestine	NM-400Y-0423	2.8 mm	2700 mm	Short	4 mm	23 G

Naming

NM-400 L-04 25

(1) (2) (3) (4)


(1) Model number

$\frac{4}{a} \frac{0}{b} \frac{0}{c}$ a (Odd: Reusable Even: Disposable)
b (0: Regular type)
c (Channel dia. 0: ϕ 2.8mm 1: ϕ 2.0mm)

(2) Sheath length

(3) Needle length (mm)

(4) Needle dia. (Gauge)

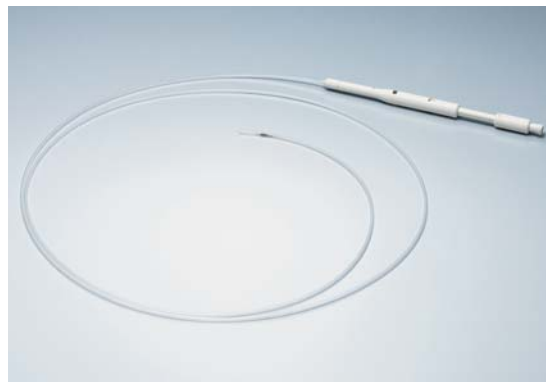

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

OLYMPUS

OLYMPUS MEDICAL SYSTEMS CORP.
Shinjuku Monolith, 3-1 Nishi-Shinjuku 2-chome, Shinjuku-ku, Tokyo 163-0914, Japan
OLYMPUS MEDICAL SYSTEMS EUROPA GMBH
Postfach 10 49 08, 20034 Hamburg / Wendenstrasse 14-18, 20097 Hamburg, Germany
OLYMPUS AMERICA INC.
3500 Corporate Parkway, P.O. Box 810, Center Valley, PA 18034-0610, U.S.A.
OLYMPUS LATIN AMERICA, INC.
5301 Blue Lagoon Drive, Suite 290 Miami, Florida 33126-2097, U.S.A.
OLYMPUS KEYMED GROUP LTD.
Keymed House, Stock Road, Southend-on-Sea, Essex SS2 5QH, England
OLYMPUS SINGAPORE PTE LTD.
481B River Valley Road #12-01/04, Valley Point Office Tower, Singapore 248373

OLYMPUS HONG KONG AND CHINA LIMITED
L43, Office Tower, Langham Place, 8 Argyle Street, Mongkok, Kowloon, Hong Kong
OLYMPUS (BEIJING) SALES & SERVICE CO., LTD.
Room 1202, NCI Tower, A12 Jianguomenwai Avenue Chaoyang District Beijing, 100022, China
OLYMPUS MOSCOW LIMITED LIABILITY COMPANY
117071, Moscow, Malaya Kaluzhskaya 19, bld. 1, fl.2, Russia
OLYMPUS AUSTRALIA PTY LTD
31 Gilby Road, Mount Waverley, VIC., 3149, Australia
OLYMPUS KOREA CO., LTD.
4F, 5F, Gyeongam Bldg., 157-27, Samsong-dong, Gangnam-ku, Seoul, Korea

www.olympus.com

E0482768